

Big Magic

Introduction

Big Magic is called “Magic” for good reason; regular use will produce truly magical results. Chronic health issues can slowly disappear, without disruption to your everyday schedule. Big Magic uses the natural Schuman frequency to induce healing. The mat covers large areas like the spine and legs, for extra convenience.

The power consumption of Big Magic is low. You can set up your wonderful Big Magic, then forget it, knowing your health will improve every night.

Pets also benefit from Big Magic. Place the mat under their sleeping pad, so they too can benefit from the naturally-healing waves.

Warning

Never use Big Magic in, near, or around water. Big Magic’s magnetic pulses also prohibit its use under the following circumstances:

- During pregnancy.

- While using a non-MRI-safe IUD.

- People or animals with organ transplants.

- People or animals with pacemakers.

- Hearing aids, medical pumps, or other electrical implants/devices.

Non-magnetic metal rods, screws, and plates are safe for PEMF use. If an implant is safe for an MRI, it is safe to use Big Magic.

Parcel Contents


Your parcel contains a single, large box.



This box has two layers. The top layer has the main unit, 5V power supply and USB power cable. The lower layer has the PEMF mat and coil tester.



Main Unit Overview

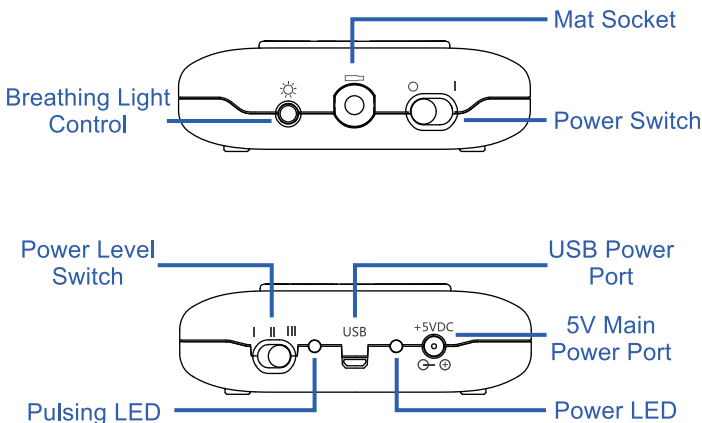

Preparing For Use

There are 4 simple steps to prepare Big Magic for use:

- a.** Plug in the power supply.
- b.** Connect and position the mat.
- c.** Set the desired power level and turn on.
- d.** Test.

Here are the steps in more details:

Plug in the power supply

There are two options to supply power to Big Magic.

- 1) Mains power, using the supplied 5v power adapter.
- 2) USB power source, such as a 5v portable power bank, 5v phone charger, or computer.



Power Adapter



**USB Power Source
(5v Phone Charger)**

You can choose the power option which best suits your needs.

Connect and position the mat

Plug the mat into Big Magic. Make sure the plug key matches the groove in the socket.




Place the PEMF mat under the area of your body where you wish to receive the PEMF treatments.



The mat is long enough to cover your entire spine. For a full night treatment, you can place it under the sheets of your bed. You may place a towel over the pad if you want greater comfort.

Set the desired power level and turn on



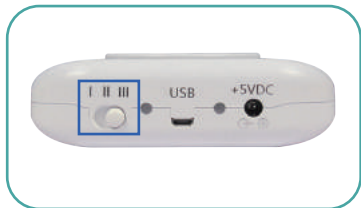

Set the power level using the 3-position switch. Level three is the one we recommend. If you are sensitive, you can select a lower power level.



Turn on Big Magic using the switch located on the right.

Test

Place the Coil Tester on the PEMF mat. The red lights should flash rapidly during operation. The coil tester should be placed in the middle of the mat like this.

Coil Tester Position

Indicator Lights

Big Magic has 3 lights:



1. Power LED

The right LED illuminates when Big Magic is receiving power.



2. Pulsing LED

The left LED flashes each time Big Magic sends an electrical pulse to the mat.



3. Breathing LED

The beautiful breathing light shows Big Magic is running.



If you place the unit on a bedside cabinet, the breathing light may distract you from sleeping. Pushing this button turns off the light.

Tips

- 1.** The PEMF mat can be folded for ease of use and transport. Convenient for business trips or family vacations.
- 2.** Clean the main unit and the mat using a damp cloth, and dry before using.
- 3.** For extra convenience, Big Magic can run continuously, with very low power consumption.